Process Challenge Device (PCD)

Quality certification

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STEAM

LOT

BT225 SCBI

Geobacillus stearothermophilus ATCC® 7953

LOT

Heat shock population / Población _____ CFU / UFC

D-value / Valor D (121 °C) (*) _____ min.

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Survival time = (log₁₀ labeled population - 2) x labeled D-value

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Kill time = (log₁₀ labeled population + 4) x labeled D-value

Z-value / Valor Z _____ °C

Parameters determined at time of manufacture according to ISO 11138-1: 2017, ISO 11138-3: 2017 and IRAM 37102: 1999 (Parts 1 and 3) standards. The values shown are reproducible only under the same conditions under which they were determined.

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ISO and USP Compliant

Terragene® and Bionova® are registered trademarks of Terragene S.A. ATCC® is a registered trademark of American Type Culture Collection.

* Determined by the fraction negative method

IT26-C Chemical Integrator

LOT

Stated Values (SV)

Temperature	Time
121 °C _____	min.
128 °C _____	min.
135 °C _____	min.

SV = stated value; Value of a critical process variable at which the indicator is designed to reach its endpoint as defined by the manufacturer.

According to ISO 11140-1:2014 standard

Lic. Adrián J. Rovetto
Technical Director

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Explanation of Symbols

Product designed for use with Steam sterilization cycles.

CE mark.

Authorized representative in the European Community.

Medical Device.

Batch number.

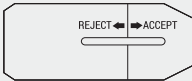
Manufacture Date.

Expiration Date.

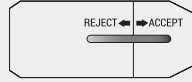
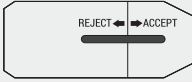
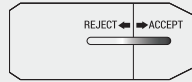
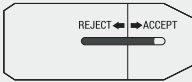
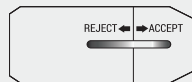
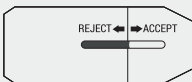
Manufacturer.

Result reference guide

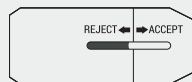
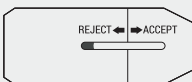
Unprocessed

Accept

Reject

EN Process Challenge Device

For a rapid and easy monitoring of Steam sterilization processes

Indications for use

United States

Terragene® Bionova® PCD (PCD225-C) provides a defined challenge resistance against the claimed cycles shown below and demonstrated resistance equivalence to the ANSI/AAMI 16 towel pack. The device provides routine monitoring and sterilizer qualification testing steam sterilization processes.

• Gravity-displacement Steam Sterilization Cycles

132 °C, 25 minutes

132 °C, 15 minutes

132 °C, 10 minutes

135 °C, 10 minutes

• Dynamic-air-removal Steam Sterilization Cycles

132 °C, 4 minutes

135 °C, 3 minutes

Outside the United States

Bionova® PCD225-C process challenge device has been designed for rapid and easy monitoring of Steam sterilization processes between 132-135 °C.

Device description

Bionova® PCD225-C Process Challenge Device (PCD) Test Pack has been designed to simulate a load to be sterilized and to pose a challenge to the sterilization process. It is used to evaluate the effective performance of the process by detecting inadequate air removal and steam penetration. It also allows release of routine loads, especially implants. Furthermore, it allows to perform the routine monitoring and periodic validation of the sterilizers (after repair, installation, relocation).

Bionova® PCD225-C consists of a disposable pre-assembled package as outlined in ANSI/AAMI ST79 which contains a BT225 Self-Contained Biological Indicator (SCBI), a PCDBI-C-RC Record Card and an IT26-C moving front chemical integrator (Type 5 according to ISO 11140-1:2014 standard) that gives instant visible indication that sterilizing conditions have been reached.

Each pack consists of a stack of porous cards holding a SCBI tube that contains a population of *Geobacillus stearothermophilus* ATCC® 7953 spores soaked on a carrier as well as growth indicator medium contained in a glass ampoule. Each SCBI has a process indicator (Type 1 according to ISO 11140-1:2014 standard) on label that changes from pink to brown when exposed to steam. The moving front chemical integrator shows ACCEPT result when sterilization conditions were reached while process indicator (Type 1 according to ISO 11140-1:2014 standard) on PCD box changes from light blue to dark grey/black when exposed to steam.

Precautions

WARNING: Do not use PCD for monitoring Ethylene Oxide, Dry Heat, Formaldehyde or any sterilization process other than Steam. Do not reuse SCBIs.

WARNING: Place one or more PCDs in sterilizing hard-to-reach areas to ensure all areas of the chamber are sterilized. Evaluate all load configurations to ensure ALL hard-to-reach areas have been identified, and place a PCD in each of those locations.

WARNING: Do not use the sterilizer until the SCBI test result is negative.

Instructions for use

1. Place the pack inside a normally loaded steam autoclave, in those areas which are considered most inaccessible for the sterilizing agent (e.g., the center of the load and areas near the door).

2. Run the sterilization cycle.

3. After the sterilization process has finished, open the sterilizer door, wait for 5 minutes and remove the test pack. **NOTE:** The color of the box may vary from the original after undergoing the sterilization cycle. This does not represent a problem regarding the operation or quality of the product.

4. Check that the process indicator printed on box has changed color from light blue to dark grey/black. Open the test pack, wait 5 minutes and remove the SCBI. Allow it to cool down to room temperature. **PRECAUTION:** Wear safety glasses and gloves when removing the SCBI from the sterilized test pack. **WARNING:** Do not crush or handle the SCBI excessively, since this might cause the glass ampoule to burst.

5. Check IT26-C moving front chemical integrator for correct

exposure. If the dark bar has reached the ACCEPT zone, this confirms that the inside of the pack has been exposed to correct sterilization conditions (refer to Result Reference Guide). Otherwise, check the sterilization process.

6. Check the process indicator printed on SCBI's label. A color change to brown confirms that the SCBI has been exposed to steam. **IMPORTANT:** This color change does not indicate that the process was sufficient to achieve sterility.

7. Identify the Bionova® BT225 SCBI by writing the sterilizer number (in case of having more than one), load number, and processing date on the label. Fill out the required information on the Record Card.

8. Press the cap to seal the tube. Crush the glass ampoule contained in the SCBI with an individual ampoule crusher or with the ampoule crusher placed within the incubator's incubation area. Then shake the tube down vigorously, with movements similar to those performed to lower the temperature in a mercury thermometer, until the medium reaches the base of the tube and soaks the spore carrier entirely. **IMPORTANT:** Use a non-sterilized SCBI as a positive control in order to ensure that correct incubation conditions were met; capability of culture medium to promote rapid growth; no alteration of spore's viability due to improper storage temperature, humidity or proximity to chemicals and also correct functioning of Terragene® Bionova® auto-reader Incubator. Both, the positive control indicator and the processed indicator, should belong to the same batch.

9. Incubate the processed SCBI and the indicator used as positive control in Terragene® Bionova® auto-reader incubators at (60 ± 2) °C to get the final fluorescence result at 7-second incubation. A 48-hour readout is optional to confirm the fluorescence result by visual color change. A positive fluorescence result (or culture medium color change after the 48-hour incubation) means failure of sterilization process. If a negative fluorescence result is obtained (or culture medium remains the original color after 48-hour incubation), it means that sterilization process was successful.

10. Record the SCBI and integrator results and adhere the self-adhesive Record Card or, alternatively, only the area containing the integrator indicator.

11. Discard the pack and the indicators immediately.

Monitoring frequency

Follow facility policies and procedures, which should specify a biological indicator monitoring frequency compliant with professional association recommended practices and/or national guidelines and standards. As the best practice and to provide optimal patient safety, Terragene® recommends that every sterilization load be monitored with an appropriate biological indicator.

Storage

Store in a dark place, at temperatures between 10-30 °C, 30-80 % relative humidity. Do not freeze. Do not store near sterilizing agents or other chemical products.

Self life

PCD225-C has an expiration date of 18 months from the date of manufacture, given by SCBI that carries, when stored at recommended conditions. Do not use after expiration date. Chemical Integrators and Process Indicators have an expiration date of 18 months when used in/on PCD. Endpoint Stability Reaction: chemical indicator endpoint shall remain unchanged for a period of not less than 6 months when stored at previously indicated conditions.

Discard

Discard the PCD after use according to your country's healthcare and safety regulations. The positive SCBIs can be autoclaved in a gravity air displacement steam sterilizer at 121 °C for 30 minutes, 132 °C for 15 minutes or 134 °C for 10 minutes; or in a dynamic air removal steam sterilizer at 132 °C for 4 minutes or 135 °C for 3 minutes.

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ATTENTION: If your chemical indicator has this logo printed on it, it means that your product is compatible with exclusive Trazanto Automatic Quality Control and Traceability System for cleaning, disinfection and sterilization processes.

This System interprets the indicator's results in an accurate and automatic manner, and stores all the process data through Bionova® Software. Use Trazanto scanner to interpret the results (see instructions for use in the manual of the device).

ATTENTION: do not write on the indicator before scanning.

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