

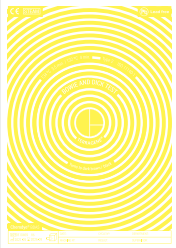
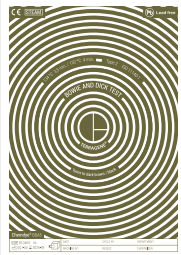
## Bowie-Dick Test Sheet


## ● Unprocessed
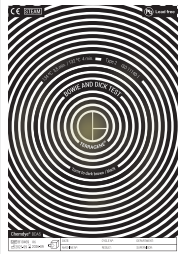
## ● Unsatisfactory
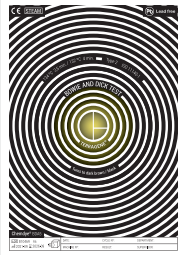
## ● Satisfactory

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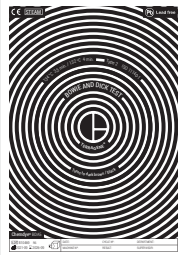Indicators 100% Lead Free

Indicators for use in  
Bowie and Dick-type test.

## Explanation of Symbols

Product designed for use in  
Steam sterilization cycles

CE mark.

Authorized representative  
in the European Community.

Manufacture Date.

Expiration Date.

Manufacturer.

NOTE: Reference colors exhibited in printed prospectus as well as those displayed in our website and soft copies of other documents, do not necessarily represent the real color as shown in the actual indicators.

## EN Bowie-Dick Test Sheet

For vacuum-assisted steam sterilizers

## Product description

Chemdye® BDA5 Bowie-Dick Test Sheets are A5 sized printed with a chemical indicator in a concentric pattern. The chemical indicator turns from yellow to dark brown/black indicating adequate air removal. An air removal failure is shown as a lighter colored area near the center of the otherwise dark brown/black Test Sheet.

## Indications for use

Chemdye® BDA5 Bowie-Dick Test Sheets are designed to be used in standard Bowie-Dick test packs assembled according to each country regulations. Test Sheets are used to test air removal efficiency at 132-134 °C in vacuum assisted steam sterilizers. The test pack should be placed horizontally on the lowest shelf in an already warmed, empty sterilizer. The sterilizer should run for 3.5 min. - 134 °C / 4 min. - 132 °C. Sterilizers which show non-uniform color development should be shut down and scheduled for service.

Chemdye® BDA5 Bowie-Dick Test Sheet is manufactured according to ISO 11140 part 1 standard.

**WARNING:** Do not use the Chemdye® BDA5 Bowie-Dick Test Sheets for air removal monitoring of gravity air displacement steam sterilization processes, Dry Heat, Ethylene Oxide, or other low temperature sterilization processes.

## Storage

Store in a dark place and at temperatures between 10-30 °C, 30-80 % Relative Humidity.

Do not freeze.

Do not store near sterilizing agents or other chemical products.

It is recommended to protect the indicator Test Sheet after use inside a plastic foil for its storage.

## Instructions for use

1. Perform the air removal test daily and after sterilizer repairs.
2. Prepare a test pack of absorbent surgical towels. The towels should be freshly laundered, but not ironed. The towels must be folded to a size not smaller than 170 mm from one side and 220 mm from the other side, and placed one above another. The height of the test pack must be 250 to 280 mm.
3. Place a single Chemdye® BDA5 Bowie-Dick Test Sheet in the center of the test pack.
4. Wrap loosely the towels with the Test Sheet and secure with self-adhesive tape.
5. Place the test pack horizontally in the front/bottom of an empty vacuum assisted steam sterilizer.
6. Run the sterilizer for 3.5 min. - 134 °C / 4 min. - 132 °C.

**NOTE:** If any serious incident occurs in relation to the device, it should be reported to Terragene S.A. and the competent authority of the State in which the user is established.

## Interpretation of results

After completing the cycle, retrieve and examine the test sheets using the Result Reference Guide.

**Satisfactory test:** This is indicated by a Test Sheet that shows uniform dark brown/black color development. Such a test indicates rapid steam penetration and adequate air removal.

**Unsatisfactory test:** This is indicated by a heterogeneous color change into dark brown/black on the Test Sheet, such as lighter areas in the center.

Unsatisfactory test should be reported immediately to the supervisor. The sterilizer should not be used until examined for possible malfunction.

## Shelf life

Chemical Indicators have an expiration date of 5 years from the date of manufacture when stored at recommended storage conditions. Both, product expiration date and batch number are specified on product packaging and in each individual card. Do not use indicators after their expiration date.

Endpoint stability: chemical indicator endpoint shall remain unchanged for a period of not less than 6 months when stored at previously indicated conditions.

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