
Infection Control Division | **Cleaning Indicators**

Lumenia L122 Cleaning Challenge Kit For Automatic Endoscope Reprocessors

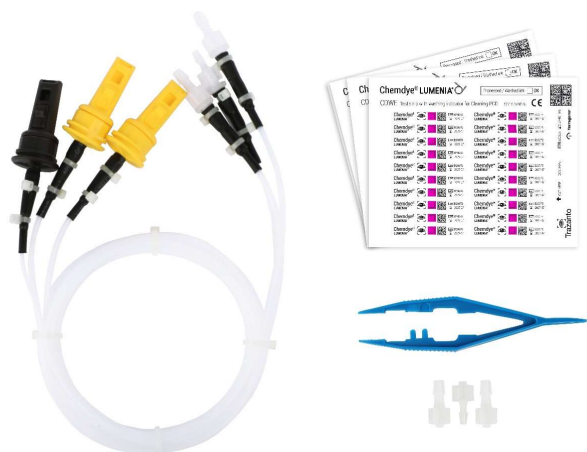Usage

Challenge Kit designed to evaluate the effectiveness of flexible endoscopes cleaning to be reprocessed in Automatic Endoscope Reprocessors.

Applicable regulation

Designed under Quality Management System standard
ISO 13485:2016/NS-EN ISO 13485:2016.
ISO 15883-4:2018, ISO/TS 15883-5:2005

Characteristics

Challenge Device

Chemdye® Lumenia L122 Kit consists of a challenge device (Chemdye® Lumenia CCDER122) and cleaning indicators (Chemdye® Lumenia CDWE) for monitoring cleaning performance in Automated Endoscope Reprocessors. The challenge device consists of three polytetrafluoroethylene (PTFE) tubes of high strength and durability, one with an inner diameter of 1 mm and two with an inner diameter of 2 mm, all of them 2 m long. Each PTFE tube has an input terminal in one of its ends for plugging it to the Automated Endoscope Reprocessor and a device in the other end that simulates an endoscope's lumen and where the Chemdye® Lumenia CDWE cleaning indicators are placed (one in the inner cavity of the device and the other one in an external slot).

Cleaning indicator

Chemdye® Lumenia CDWE cleaning indicator was especially designed for monitoring these processes. It consists of a test soil, developed in accordance with the requirements described in ISO/TS 15883-5:2005 standard, arranged on a self-adhesive thermolabile synthetic support. The formulation of the cleaning indicator allows to detect the factors that affect the result of the cleaning process.

Size of the indicator strip: 38 mm x 6 mm.

Initial condition: fuchsia color.

Final condition: total elimination of the fuchsia organic substance.
The indicator's ink is not toxic.

NOTE: If unfavorable results are obtained, check quality/quantity of chemicals, water pressure, temperature, cleaning time and general condition of the washing machine.

Environmental conditions during manufacture

T = 15-30 °C, R. H. = 30-80 %.

Storage conditions


T = 10-30 °C, R. H. = 30-80 %. Store in dark place.

Transport conditions

Storage conditions should be strictly followed.

Infection Control Division | **Cleaning Indicators**

Lumenia L122 Cleaning Challenge Kit For Automatic Endoscope Reprocessors

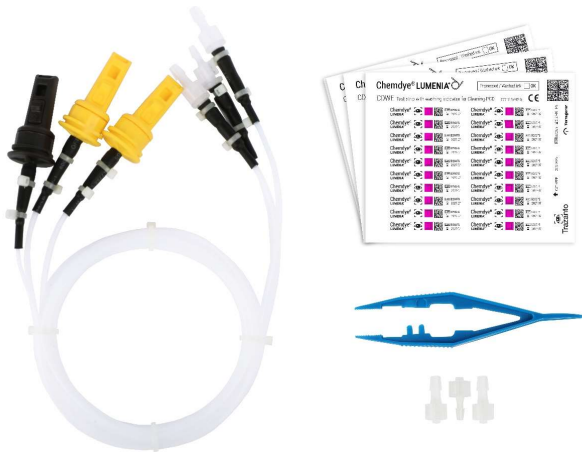Transport in closed and strengthened boxes in order to avoid damages.

The transport of this product does not represent a risk for health.

Shelf life

Cleaning indicator: 5 years.

Packing

1 Challenge Device + 18 sheets of cleaning indicators (360 indicators) + 3 input terminals + 1 plastic tweezer + Instructions for Use.

Packing information: product code and description, intended use, regulation, storage conditions, manufacturer information and data on label.

Labelling

On product's pack: product code, kit batch number and its components' batch numbers, presentation, manufacture and expiration date, barcode and datamatrix code.

Possible target markets

Healthcare and Medical Industries

Other important information

Read product's instructions for use thoroughly before use.

Precautions

Store the product within its original package, away from light, disinfectants, gases and humidity.

Do not re-use the indicator.

Do not touch indicator ink with your hands.

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